


The NPS Anti-Static Mobile Kit

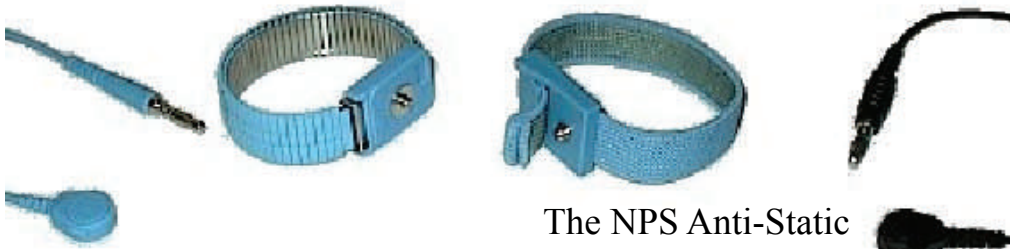


The NPS Anti-Static Mobile Kit improves data acquisition results and protects your equipment.

The NPS Anti-Static Kit is easy to use, cost effective and of the highest quality.

Noticeable features:


- Easy setup design
- Metal wrist band for the best connection
- Adjustable wrist band for those with really small or large hands and wrists
- Long cords a more comfortable reach
- Handy storing snaps for disconnected cables

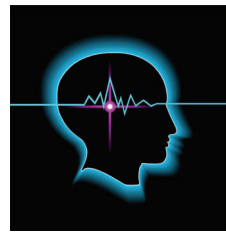

NeuroPulse-Systems LLC


324 E. Railroad Ave., Suite 400
Ft. Morgan, CO 80701
(970) 867-2527 / Fax: (775) 218-3241

www.np-systems.com

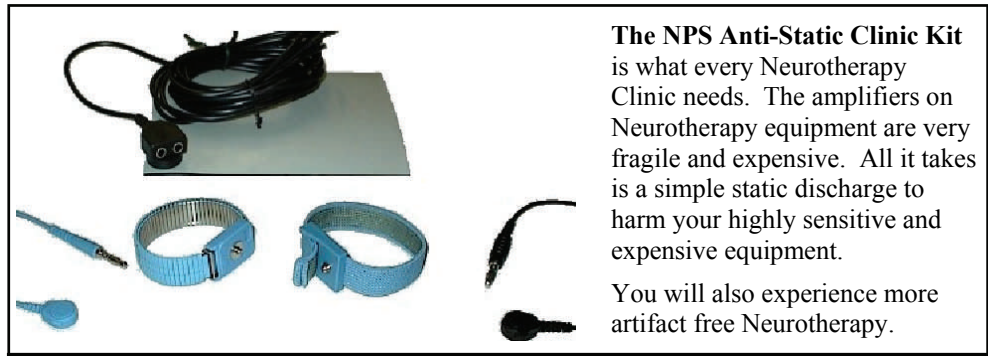


The NPS Anti-Static Clinic Kit is what every Neurotherapy Clinic needs. The amplifiers on Neurotherapy equipment are very fragile and expensive. All it takes is a simple static discharge to harm your highly sensitive and expensive equipment.

You will also experience more artifact free Neurotherapy.

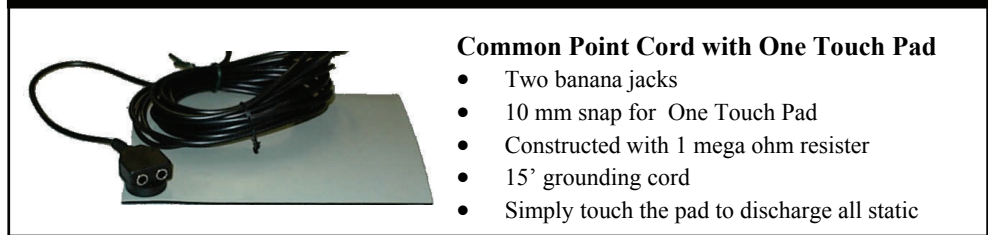


Common Point Cord with One Touch Pad

- Two banana jacks
- 10 mm snap for One Touch Pad
- Constructed with 1 mega ohm resistor
- 15' grounding cord
- Simply touch the pad to discharge all static

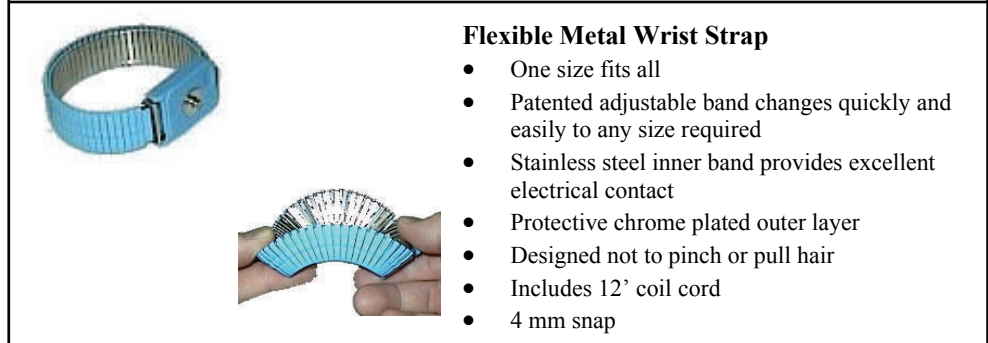


Flexible Metal Wrist Strap

- One size fits all
- Patented adjustable band changes quickly and easily to any size required
- Stainless steel inner band provides excellent electrical contact
- Protective chrome plated outer layer
- Designed not to pinch or pull hair
- Includes 12' coil cord
- 4 mm snap

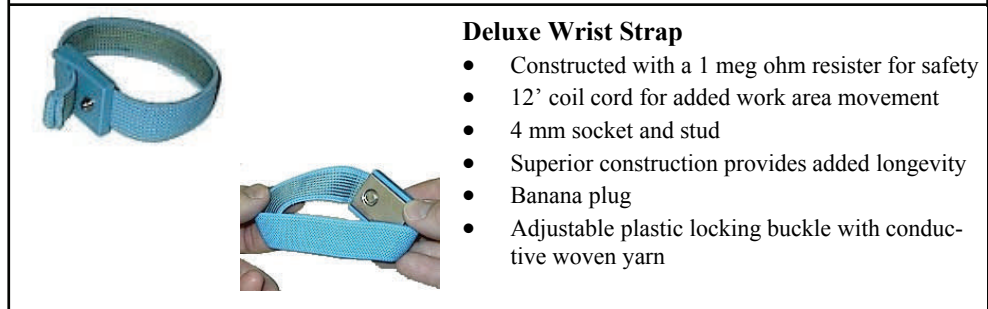

Deluxe Wrist Strap

- Constructed with a 1 meg ohm resistor for safety
- 12' coil cord for added work area movement
- 4 mm socket and stud
- Superior construction provides added longevity
- Banana plug
- Adjustable plastic locking buckle with conductive woven yarn

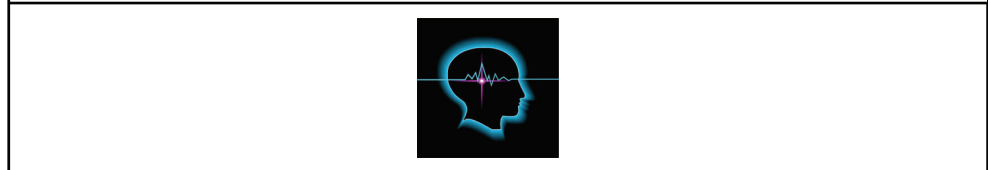

The NPS Anti-Static Clinic Kit does not protect against all types of static or electrical discharge. The operator takes full responsibility to properly check for ground, install the base station and have themselves and anyone else working with the equipment properly grounded before touching the equipment. Good connections must be kept continuously for proper function. NeuroPulse-Systems LLC will not be held liable for any equipment damage, due to the fact the operator must use proper grounding techniques with or without Anti-Static equipment. The NPS Anti-Static Clinic Kit simply reduces the chances of static discharge that could result in equipment damage. It is also from our experience that Anti-Static equipment improves data results when the subject is properly grounded.

Copyright © 2004 NeuroPulse-Systems LLC. All Rights Reserved. NP24 is a trademark of NeuroPulse-Systems LLC. Features and specifications are subject to change without notice. Weights and measurements are approximate.

NeuroPulse-Systems LLC makes no claims as to potential medical or therapeutic benefits.